


Switch Plate - GoPro Compatible

USER QUICK GUIDE

Introduction

This mounting plate is ideal for individuals looking for a cost-effective way of mounting an assistive switch and hoping to utilize readily available GoPro mounting accessories.

Features

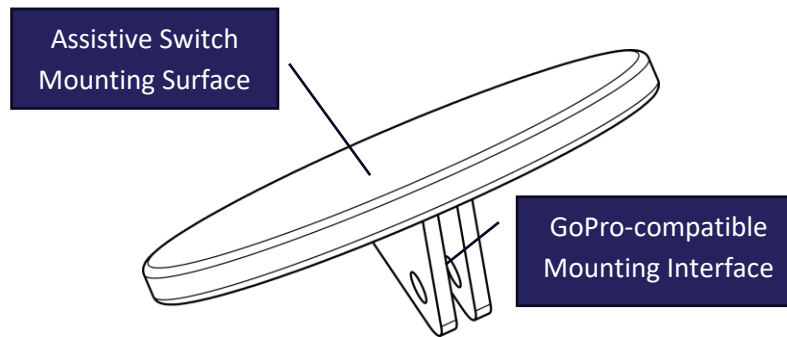

Specifications

Diameter	Ø89 mm (Ø3.5 inch)
Mass	25 g (Switch Plate Only)

Compatibility

Assistive Switches

The Switch Plate – GoPro Compatible is compatible with Assistive Switches that can be mounted to a flat surface. This includes commercial options like AbleNet Jelly Bean and DIY OpenAT options like the Interact Switch and MMC 60 Switch.

Mounting

The Switch Plate – GoPro Compatible is compatible with mounting accessories that have a three-fingered mounting adapter.

Switch Plate - GoPro Compatible

USER QUICK GUIDE


Usage

Connecting to GoPro Mounting Accessory

The two-finger mounting interface on the back of the Switch Plate is inserted into the three fingers of the GoPro Mounting Finger Adapter and secured with the GoPro Thumb Screw.

Attaching Assistive Switch

Attach the Assistive Switch to the Mounting Surface of the Switch Plate using hook and loop fastener (e.g., VELCRO) or double-sided adhesive.

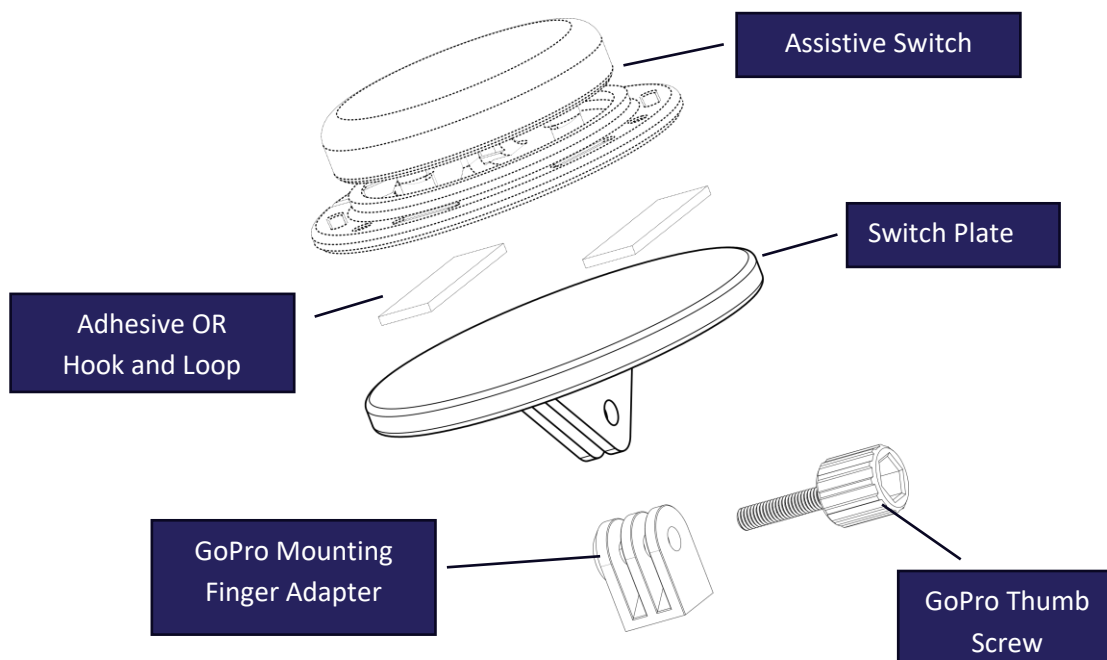

Cleaning

The Switch Plate can be cleaned with a damp cloth. The Switch Plate can also be cleaned by scrubbing with warm water and dish soap. Do not expose to hot water or a dishwasher. Do not use harsh cleaners.